


Sample Name: Poly(dimethyl siloxane -b- acrylic acid)

Sample #: P6442-DMSAA

Structure:

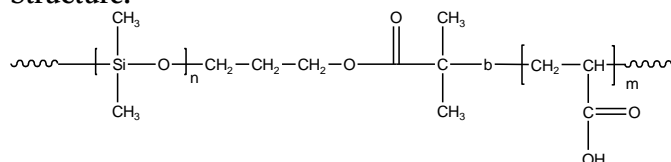


Composition:

| $M_n \times 10^3$ PDMS-b-PAA | PDI |
|---------------------------------|------|
| 8-b-11.3 | 1.80 |

Synthesis Procedure:

Poly(dimethyl siloxane-b-acrylic acid) is prepared by living anionic polymerization of hexamethylcyclotrisiloxane followed by controlled radical polymerization of protected acrylic acid monomer. The reaction scheme is shown below:

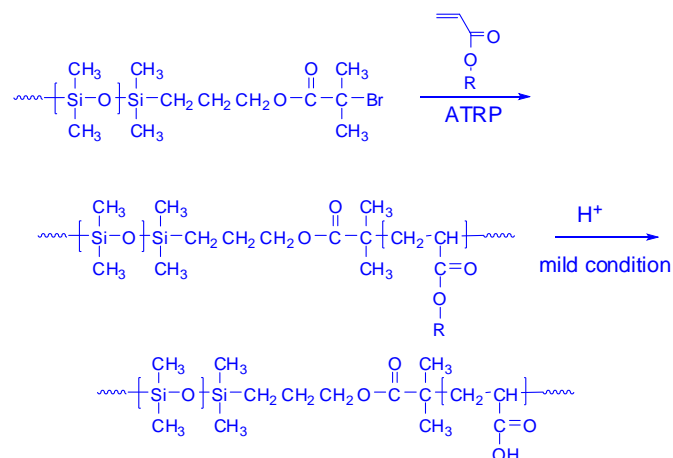

Characterization:


An aliquot of the Poly(dimethyl siloxane) block was terminated before controlled radical polymerization of the acrylate block and analyzed by size exclusion chromatography (SEC) to obtain the molecular weight and polydispersity index (PDI). The final block copolymer composition was calculated from $^1\text{H-NMR}$ spectroscopy by comparing the peak area of the siloxane protons at about 0 ppm with the peak area of the acrylic acid protons at about 1-2.5 ppm. The composition is also calculated according to FTIR showed as followed graph. The results from NMR and FTIR are compatible.

Note: The protected polyacrylic acid may cause the SEC profile broadening. We still claim the M_w/M_n as the apparent molecular weight. Real M_w/M_n should be narrower


Solubility:

Poly(dimethyl siloxane -b- acrylic acid) is soluble in THF, chloroform and dichloromethane. It is not reprecipitated well in methanol or hexanes because of its amphiphilic character.

$^1\text{H-NMR}$ Spectrum of the block copolymer:

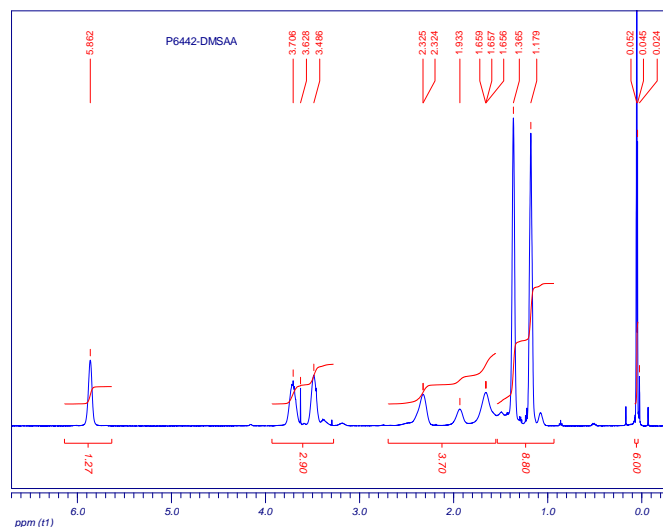

SEC profile of the block copolymer:

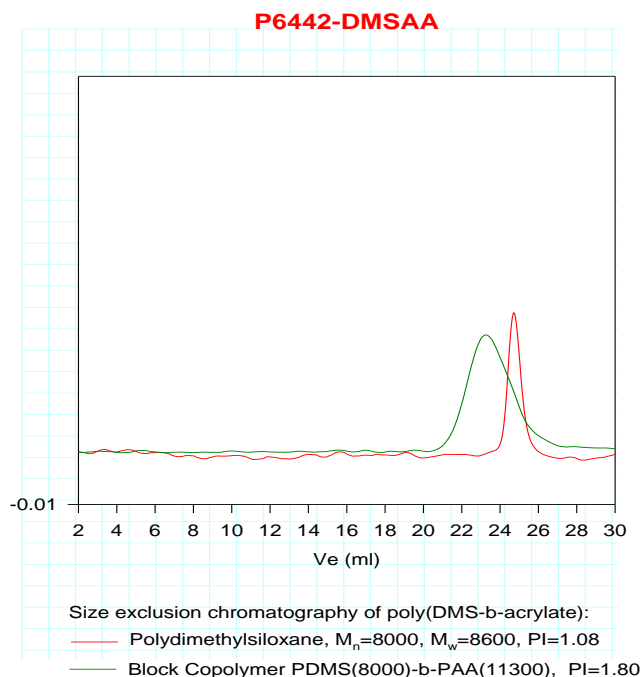


Thermal analysis of the sample# P6442-DMSAA

Thermal analysis of the samples was carried out on a TA Q100 differential scanning calorimeter at a heating rate of 10°C/min. The midpoint of the slope change of the heat flow plot of the second heating scan was considered as the glass transition temperature (T_g).

Melting and crystallization curve for the sample

The melting temperature (T_m) was taken as the maximum of the endothermic peak where as the crystallization temperature (T_c) was considered as the minimum of the exothermic peak.


Thermogram for AA block:

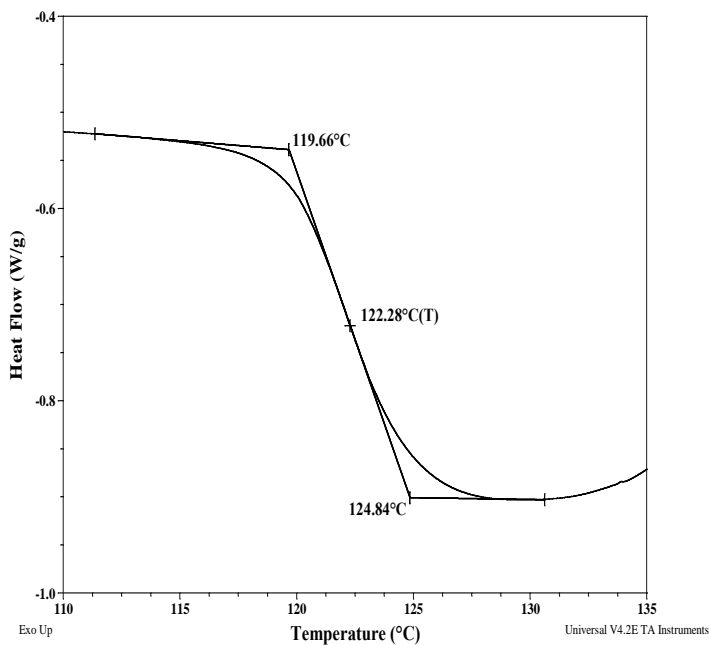

Thermal analysis results at a glance

| Sample | T_m (°C) | T_c (°C) | T_g (°C) |
|--------|------------|------------|------------|
| DMS | -44 | -71 | -127 (Lit) |
| AA | - | - | 122 |

Melting curve for DMS block:

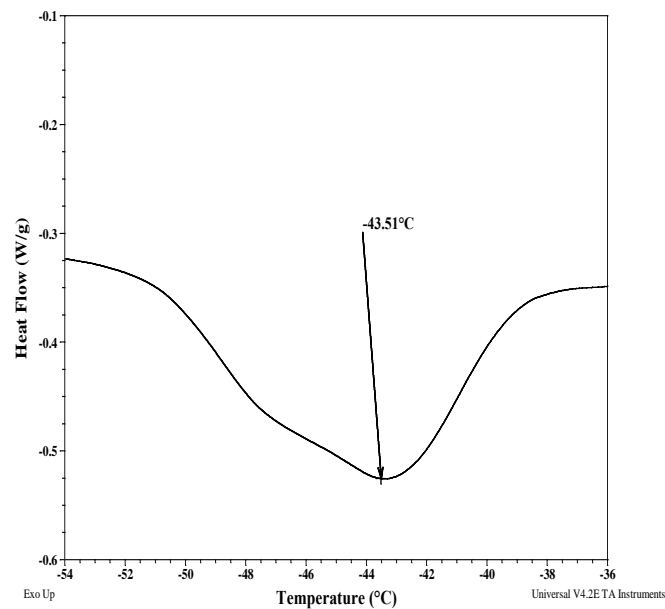

Crystallization curve for DMS:

